


Evident Fluoview 4000 (Irchel) - Shut Down

How to shut down the Evident Fluoview 4000 (Irchel) located at the Center for Microscopy and Image analysis, UZH, Room Y42-H-93.

Written By: Urs Ziegler

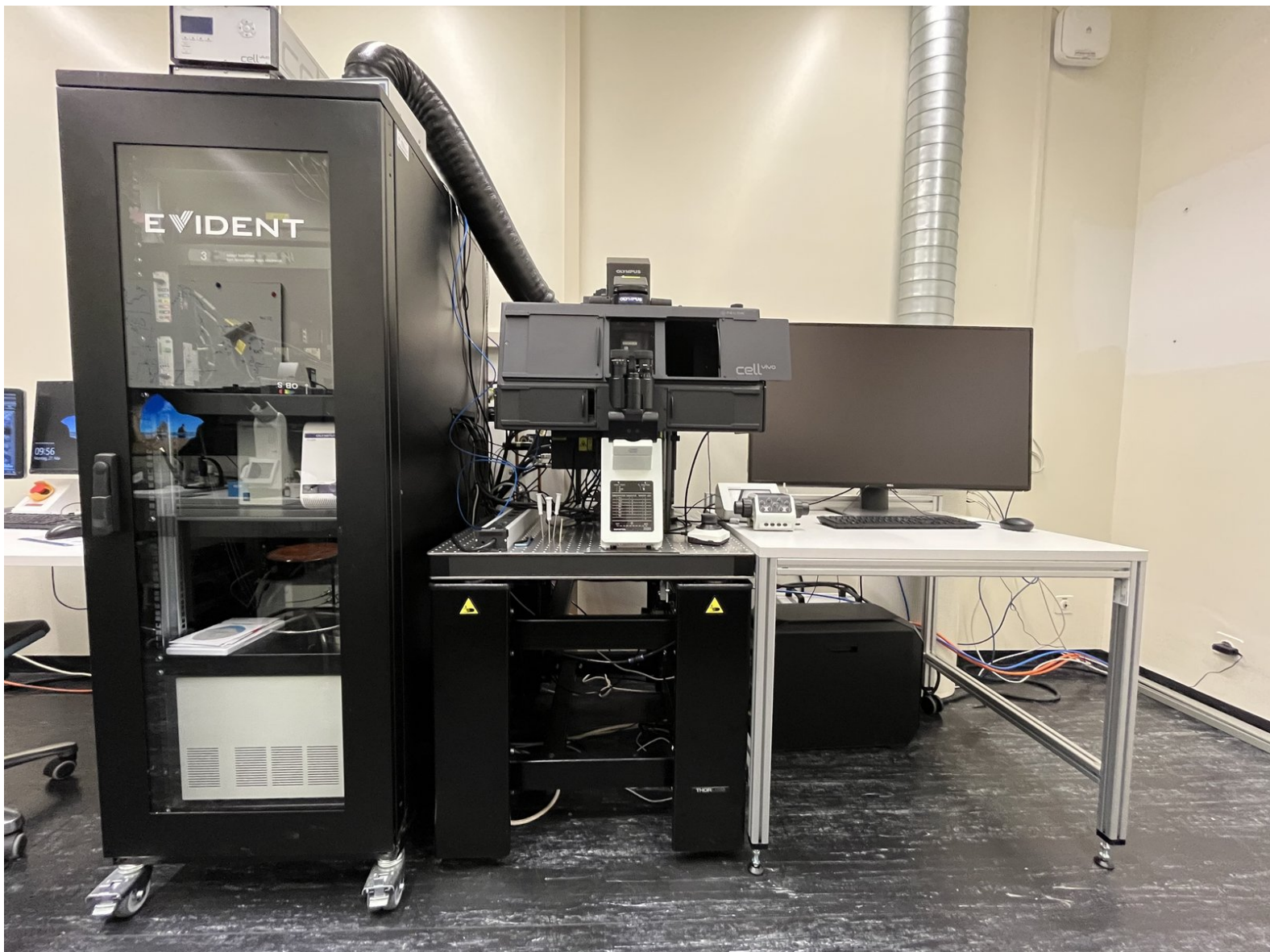


Step 1 — Save your data and close CellSens




- Confirm that you have saved your data in the appropriate location.
- Remove your sample.
- ⚠ Clean the objectives using the available Kimtech wipes and ethanol.
- Close Cellsense FV.
- If you have used the environmental control please empty the water bottle.
- ☑ Check in the booking system if there is another booking withing the next 2h:
 - If yes, you just need to log off.
 - If you are the last user of the day follow the next steps.

Step 2 — Shut down the computer, lasers and scanner

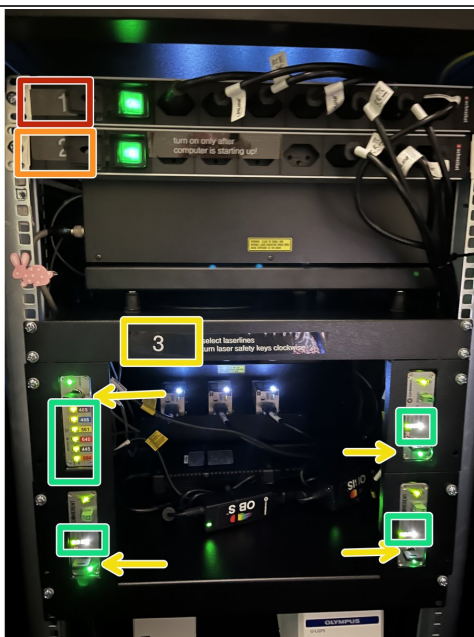

- Shut down the computer.
- Switch off the laser power strip [3].
 - ⓘ Safety interlock keys can be turned counter clock wise. No need to turn off the laser lines with the individual switches.
- Switch off the scanner unit [2].
- Switch off the main power [1]

Step 3 — Shut off environmental control

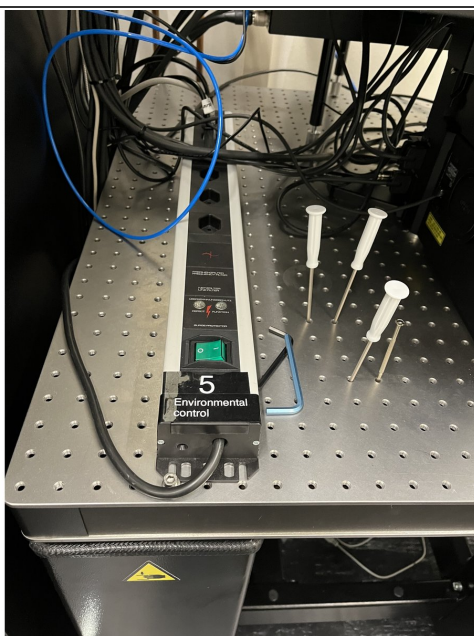

- Switch off the environmental control with the power strip [5].