


Visitron Spinning Disk - Shut down

Written By: Jana Döhner




Step 1 — Cleaning and closing software

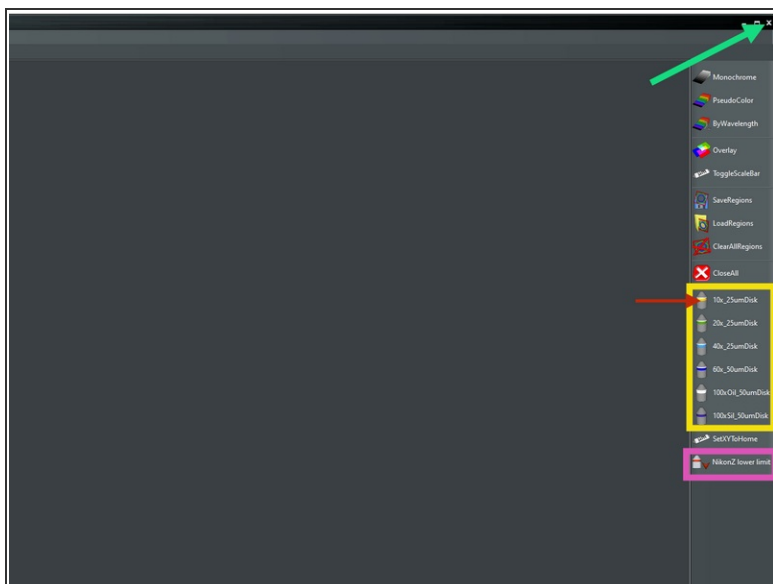



- **Lower objectives** to the lowest limit by pressing "NikonZ Lower Limit".
- Toggle between the objectives to **Clean** the used immersion objectives using the available **KIMTECH wipes** and **100% EtOH**.
- Finally select the lowest magnification - here **10x**.
- **Close** the software by crossing it out

Step 2 — Stage Insert

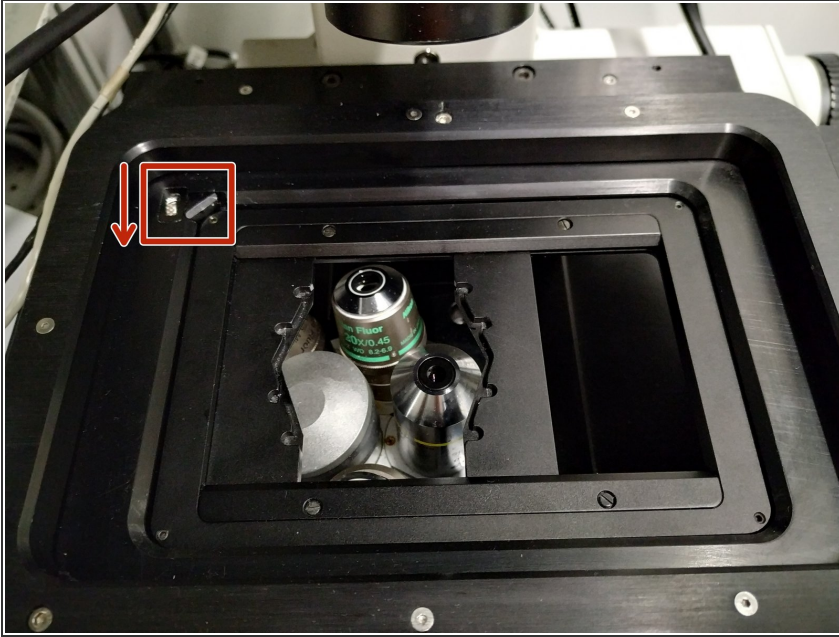

- **Remove** the used **stage insert** by releasing the screw (arrow direction indicates how to loosen the screw) and **stow away** accordingly

Step 3 — Log-Off or Shut Down

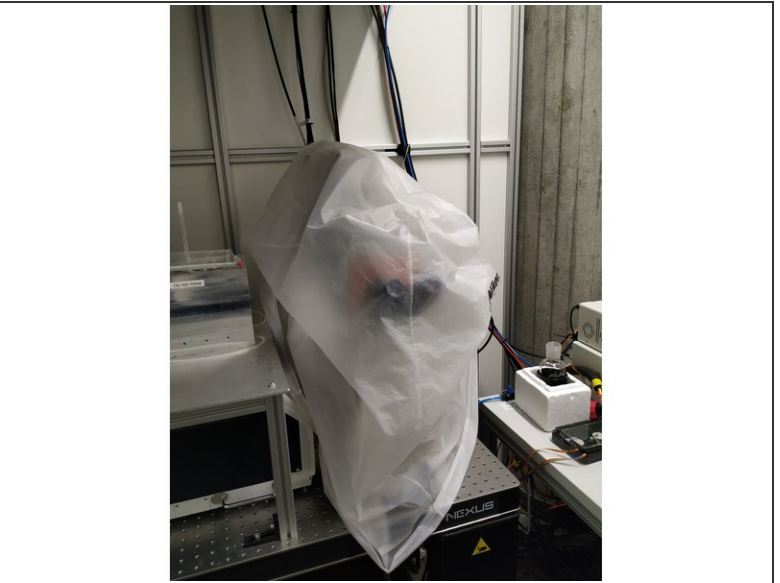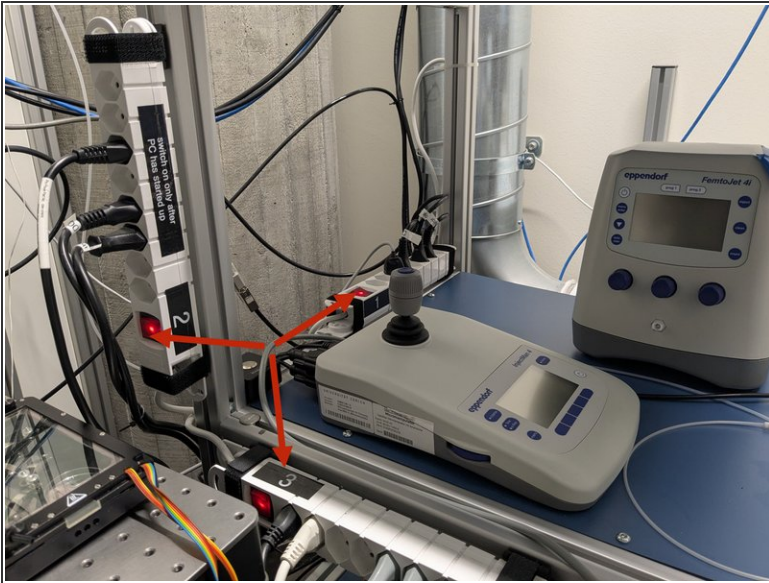

- ① Check in the **reservation system** if there is another booking shortly after you.
- If **"YES"** just sign-out from the PC.
- If **"NO"** shut-down the PC.
- After the PC has been shut-down completely, switch OFF all power supplies ([3], [2] and [1]).
- **Cover** the microscope.