Moku:Go Quickstart

A general overview on how to get the Moku:Go up and running

Written By: Angelo Regec





Step 1 — Power and Setup



- The Moku:Go can be powered by either the USB-C cable (Right) or the magnetic power cable (Left)
- If you are using the USB-C cable or both power options, go to step 3. If you are using only the magnetic power cable, see step 2.

Step 2 — Connecting the Moku:Go Wirelessly

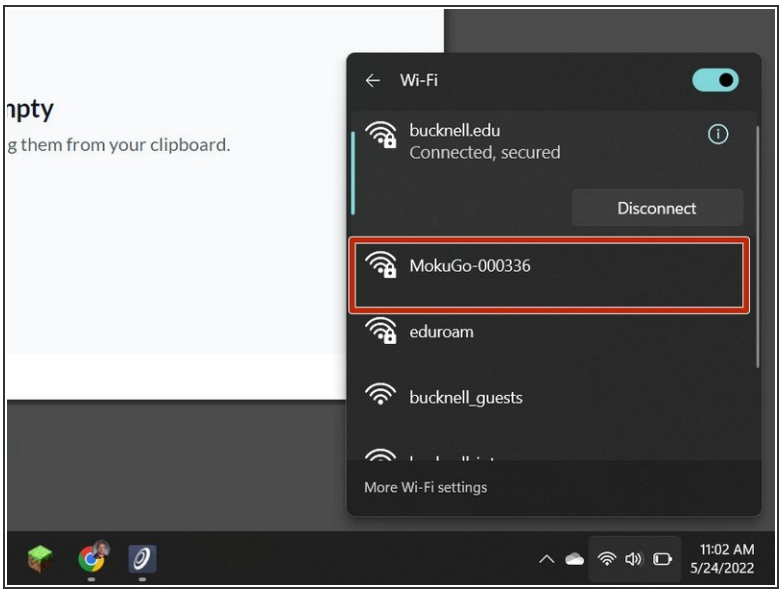


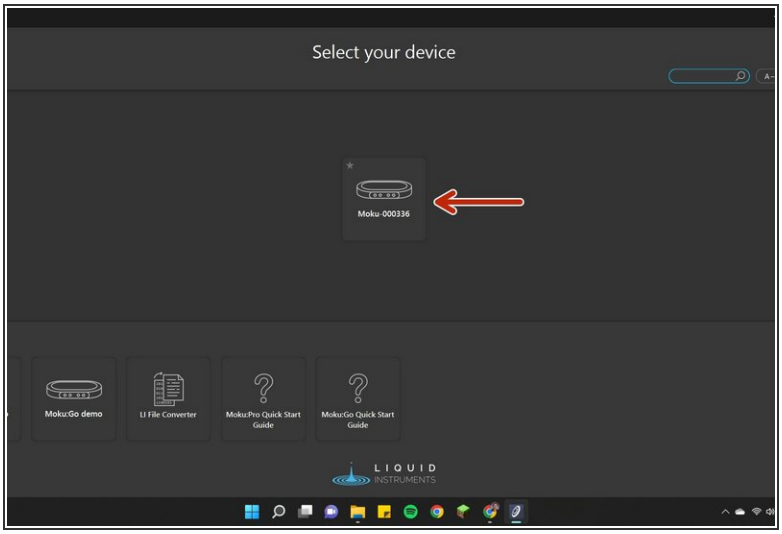


- First, open your network settings. Find the network titled "MokuGo-XXXXXX". To connect, enter the password, "password"
- Once connected, open the Moku software on your device. You should see your device shown in the Moku software as shown
- After selecting your device, the instruments should be displayed as shown.

⚠ While connected to Moku's network, you cannot access the internet. Also, you must stay connected to the Moku network for the duration of use.

Step 3 — Connecting the Moku:Go with USB-C

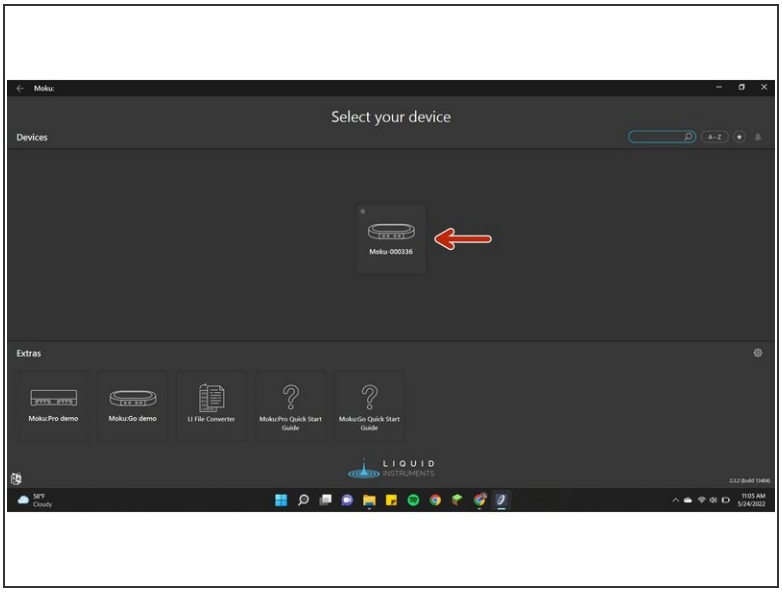

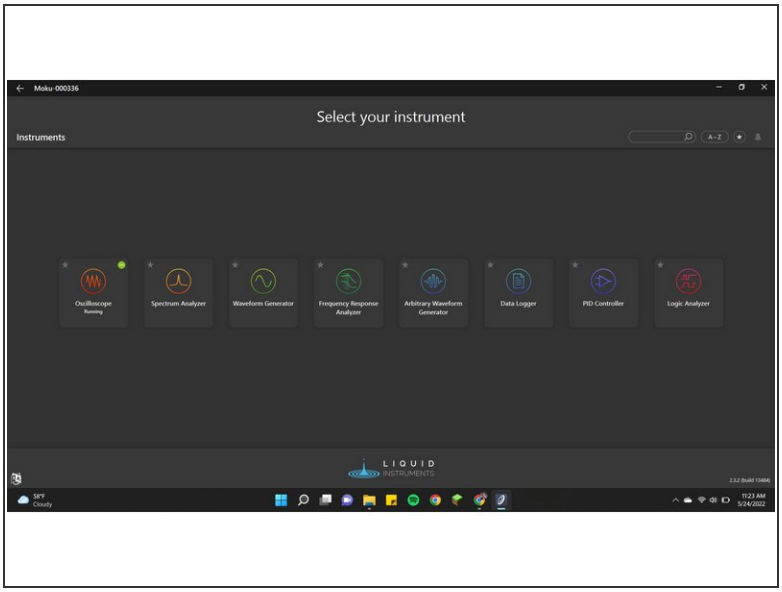

- Wait for the light on your Moku to turn green. Once it does, open up the Moku software on your device. Pictured is what you should see.
- After selecting your device, the instruments should be displayed as shown.